


## FOREHEAD THERMOMETER FLU PACK

These CE registered strips are made from flexible polyester and printed with six temperature indicating sections from 35°C to 40°C in 1°C increments. Each temperature event shows a green colour when the indicated temperature is reached; it has been designed only to be used on the forehead and has been adjusted to reflect core body temperature. When placed on the forehead, the thermometer is an accurate, fast and easy to use non-invasive method for monitoring temperature trends. Simply hold the thermometer in the centre of the forehead and you will see the green square move along the scale. After 15 to 20 seconds the temperature will settle and you can read the thermometer by looking for the green square. The forehead thermometer is reusable, unbreakable, latex and mercury free making it ideal to use in a school environment. It can be used as a simple and effective way to show changes in temperature before and after exercise. The forehead thermometer is supplied with a useful instruction leaflet including the differences in symptoms of cold vs. flu, all in a plastic wallet.

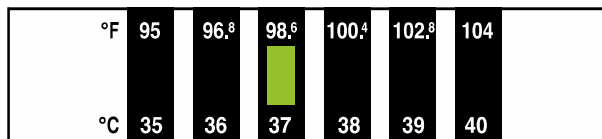

Part No.	°C	Description	Specification	Quantity
L06THSHZ035D01PK	35°C, 36°C, 37°C, 38°C, 39°C, 40°C,	Forehead Cold/Flu Pack	Thermometer, instructions and plastic wallet	Pack of 10

## Special Requests

If you have an idea for use in education that involves colour change then please contact us by going to [www.lcrhallcrest.com/thermosmart](http://www.lcrhallcrest.com/thermosmart) and clicking on the Special Request area